

## ADVANCED CATHETERISATION TRAINER - COMPLETE SET

**60150** COMPLETE SET


Simulated Patient

Our Advanced Catheterisation Trainer facilitates the teaching of urethral and suprapubic catheterisation, as well as how to demonstrate self catheterisation to patients (when using the optional Self Catheterisation Stand). It has interchangeable accurate male and female anatomy, and both sets of genitalia are realistically soft and pliable. The feel of the catheter passing along the urethra into the bladder corresponds closely to real life.

Management of suprapubic catheters can also be taught, using our optional pre-ported Suprapubic Catheterisation Management Bungs.

### SKILLS

- Correct handling of male and female anatomy
- Aseptic catheterisation technique
- Catheter placement: accepts catheters from 14-16 French (male) and 12-16 French (female)
- Fluid management
- Withdrawal of catheter
- Intermittent Self Catheterisation (ISC) when used with the optional Self Catheterisation Stand (**60167**)
- Suprapubic catheter insertion using Suprapubic Bung
- Suprapubic catheter changing using optional pre-ported Suprapubic Catheterisation Management Bung (**60168**)

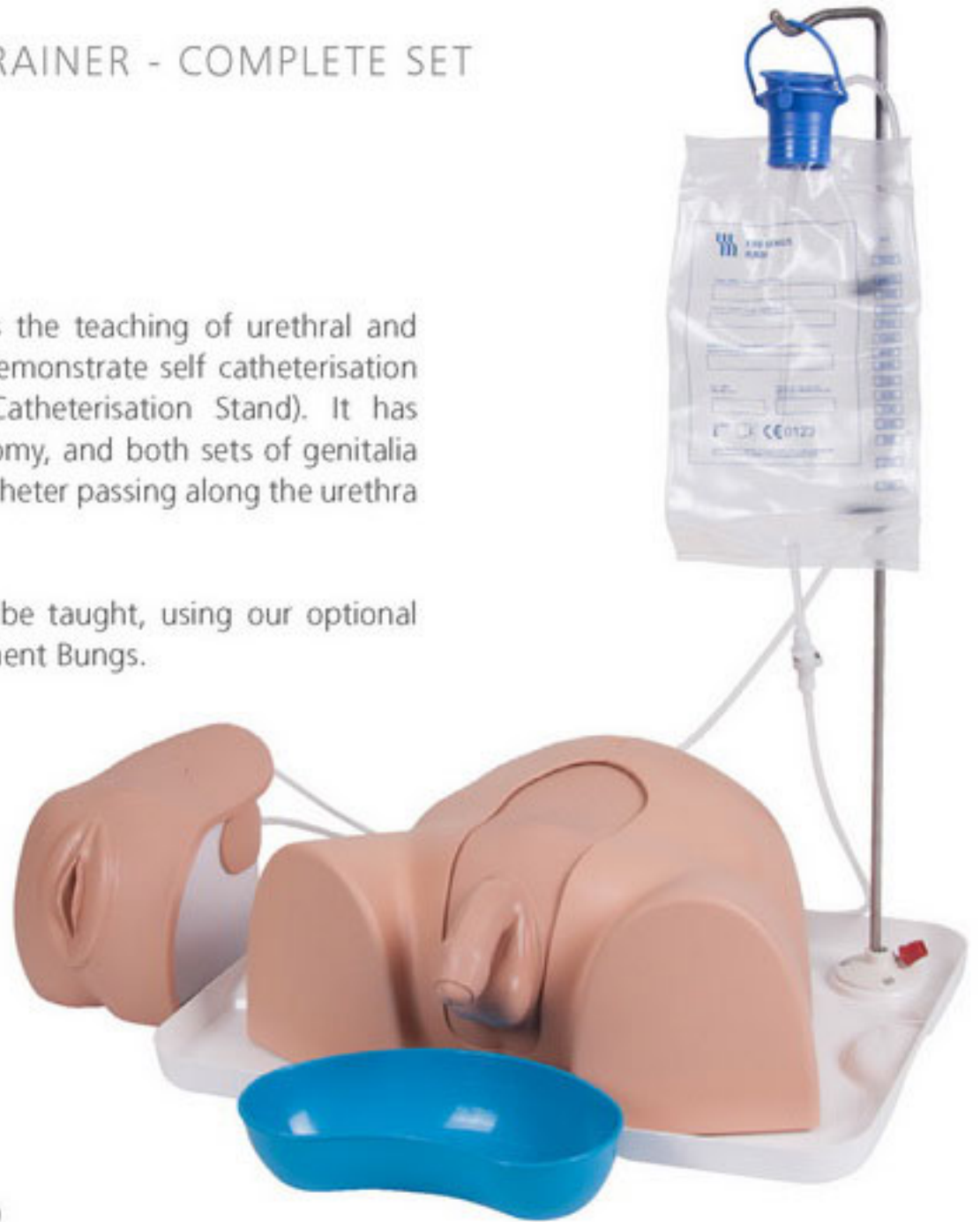

### QUALITIES

#### ANATOMY

- Male anatomy: realistic meatus, flaccid penis and replaceable foreskin
- Female anatomy: partable labia with realistic urethral meatus and vaginal opening

#### REALISM

- Supple urethra and resistant sphincter providing realistic response

#### KEY FEATURES

- 1500 ml fluid bag on stand provides pressurised fluid flow
- Syringe supplied with water-based lubricant to simulate proprietary local anaesthetic gel
- Non-drip valve
- Reusable double-sleeve catheter packaging is supplied for teaching aseptic catheter handling

#### VERSATILITY


- Easily interchangeable male and female catheterisation units

#### SAFETY

- Product is latex free (catheters provided contain latex)

#### CLEANING & CARE

- Skin surface is washable using soap and water

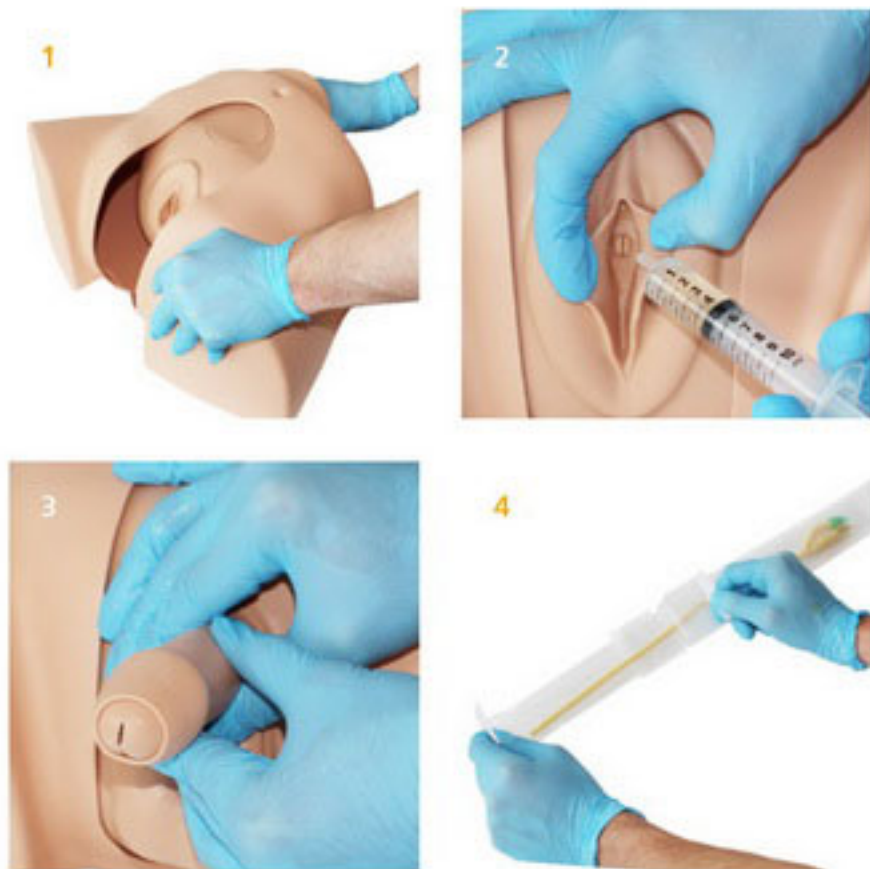

- 1 Easily interchangeable units
- 2 Accurate anatomical modelling
- 3 Penis has replaceable foreskin
- 4 Aseptic catheterisation technique

## ADVANCED CATHETERISATION TRAINER - MALE & FEMALE

The Advanced Catheterisation Trainers are available separately. Both are supplied with a suprapubic site and bung.



## BASIC CATHETERISATION TRAINERS - VALUE SET

### 60166 VALUE SET

A cost effective Value Set containing basic versions of each of the Male and Female Catheterisation Trainers. These do not contain the suprapubic bungs.



### OPTIONAL

## SELF CATHETERISATION STAND

### 60167

Use with Catheterisation Trainers: **60150, 60151, 60155, 60166** for teaching intermittent self catheterisation.



## SUPRAPUBIC CATHETERISATION MANAGEMENT BUNG

### 60168

This Pre-Ported Bung allows for the management of catheters after surgical procedure has taken place.

Use with Catheterisation Trainers: **60150, 60151, 60155.**



## CONTAINS

	60150	60151	60155	60166
• 60152 Male Catheterisation Unit	✓	✓	✗	✓
• 60154 Catheterisation Foreskin (x3)	✓	✓	✗	✓
• 60156 Female Catheterisation Unit	✓	✗	✓	✓
• 60158 Suprapubic bung (x2)	✓	✓	✓	✗
• 60161 Standard bung	✓ (x2)	✓	✓	✓ (x2)
• 60159 Male Catheterisation Sphincter Kit	✓	✓	✗	✓
• 60162 Female Catheterisation Sphincter Kit	✓	✗	✓	✓
• 60160 Aseptic Catheterisation Sleeves	✓ (x2)	✓ (x2)	✓ (x2)	✓ (x4)
• 50183 Mock LA Gel with Syringe	✓	✓	✓	✓ (x2)
• 50185 Antifungal Agent (Milton Sterilising Fluid)	✓	✓	✓	✓ (x2)
• Drip tray with bag stand	✓	✓	✓	✓ (x2)
• 1500 ml fluid bag & plastic kidney dish	✓	✓	✓	✓ (x2)
• 14 French foley catheter (male)	✓	✓	✗	✓
• 12 French foley catheter (female)	✓	✗	✓	✓
• Pelvic shell	✓	✓	✓	✓ (x2)
• 2 Litre jug	✓	✓	✓	✓ (x2)
• Carry case	✓	✓	✓	✓ (x2)