

# Injection and Blood Sampling Practice Model Type I

with storage case

## Outline

A series of techniques such as confirmation of injection site, insertion hypodermic needle and injection of medicine are possible with this model. The pressurized blood vessel can be palpated, and needle insertion provides a feeling similar to that offered by a real human arm. Blood vessel tube is durable enough to be used in numerous practice injections.



Needle injection has same feel as with a human arm.



Replaceable outer skin and blood vessel tube.




## Features

- Made of a special silicone rubber, this lifelike model limb features synthetic skin that is remarkably similar to human skin.
- The sensation of injecting a hypodermic needle into this model's simulated skin and blood vessels is identical to that encountered when performing an injection on an actual patient. Practice of both medicine injection and blood sampling using this model provides trainees with a remarkably realistic experience.
- Resilient synthetic skin conceals injection marks well, even after numerous practice injections. Injection points are highly resistant to tears and other damage.
- The artificial blood features the same viscosity and color as those of actual human blood. Since it is entirely water soluble, spills are cleaned easily without the worry of staining.
- The model's lifelike circulation operates through a simple manual process. This rugged, serviceable model incorporates extremely reliable mechanisms for trouble-free service over the long term.
- This simply designed model enables easy exchange of both the synthetic skin and veins.

## Components

- ① Arm skin
- ② Arm musculature
- ③ Blood vessel tube
- ④ Stand
- ⑤ Irrigator
- ⑥ Simulated blood
- ⑦ Elbow rest

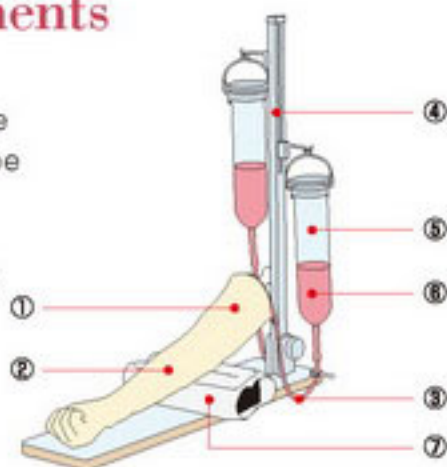

### Accessory Set LM-028A

1pc of arm skin, 5pcs of blood vessel tube, 500ml of simulated blood

### Sold Separately

LM-028H	Plastic container	2
---------	-------------------	---



## Specifications

Arm	Approx.59(L) × 9(W) × 9(H)cm	Approx.2kg
Stand	Approx.60(L) × 16(W) × 67(H)cm	Approx.1.5kg
Storage case	Approx.67(L) × 30(W) × 21(H)cm	Approx.4kg

## Components

Arm skin	1
Arm musculature	1
Blood vessel tube	6
Base and stand	1
Irrigator (500ml)	2
Simulated blood (500ml)	1
Elbow rest	1
Pinchcock	2
Tube connector	2
Baby powder	1
Transfusion set	1
Storage case	1

## Spares

LM-028B	Arm skin	1
LM-028C	Blood vessel tube	10
LM-028D	Simulated blood (500ml)	1
LM028E	Tube connector	2
LM-028S-T	Arm Skin (Thick type)	1